



## **Biocard**<sup>®</sup> MALARIA Pf RAPID TEST

Biocard<sup>®</sup> Malaria Pf Rapid test kit is a immunochromatography, qualitative analysis test kit for detection of *Plasmodium falciparum* antigen in whole blood. It is coated with target anti-HRP-II (Histidine Rich Protein) for *Plasmodium falciparum* detection



Technique: Rapid test based on lateral flow



Accuracy: Meets WHO recommended detection criteria



Detection: *Plasmodium falciparum* (Pf)



Results: 20 minutes



Longer shelf life: 24 months



Storage: 4°C to 30°C



Convenient pack size: 30 tests



Sensitivity: 100%



Specificity: 100%


## Easy to Use

- Simple test using only one step reagent
- Ease of use, requires minimal training of personnel
- Requires only 5 µL of venous or capillary whole blood
- No high cost equipment required for testing
- No need for skilled microscopist


## Easy to Interpret

- Meets WHO acceptance criteria
- Excellent correlation with microscopic results
- Differentiated results allow diagnosis of potentially lethal *P. falciparum* infection
- Easy to read, qualitative results
- Higher Specificity and better Sensitivity

## Test Procedure

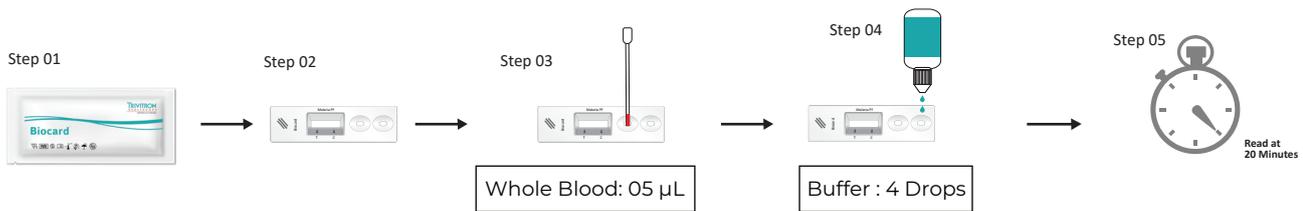

## Result Interpretation

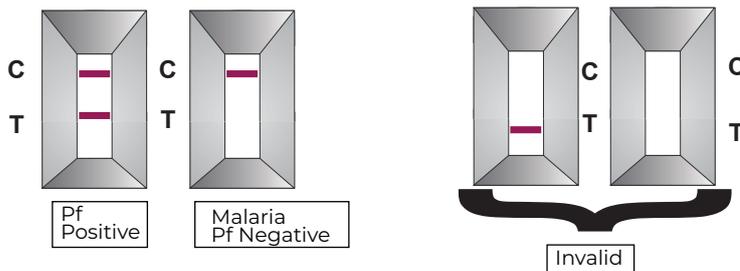

## Ordering Information

Catalog	Product	Pack Size	Specimen
30201434	Biocard® Malaria Pf Ag Rapid Test	30 Tests	Whole Blood

\*Please note: Product specifications are subject to change without prior notice owing to product modifications, improvements / up-gradation.



corporate@trivitron.com | www.trivitron.com

Customer Care Helpline: +91 98400 80008



**TRIVITRON**  
HEALTHCARE  
speaking your language